
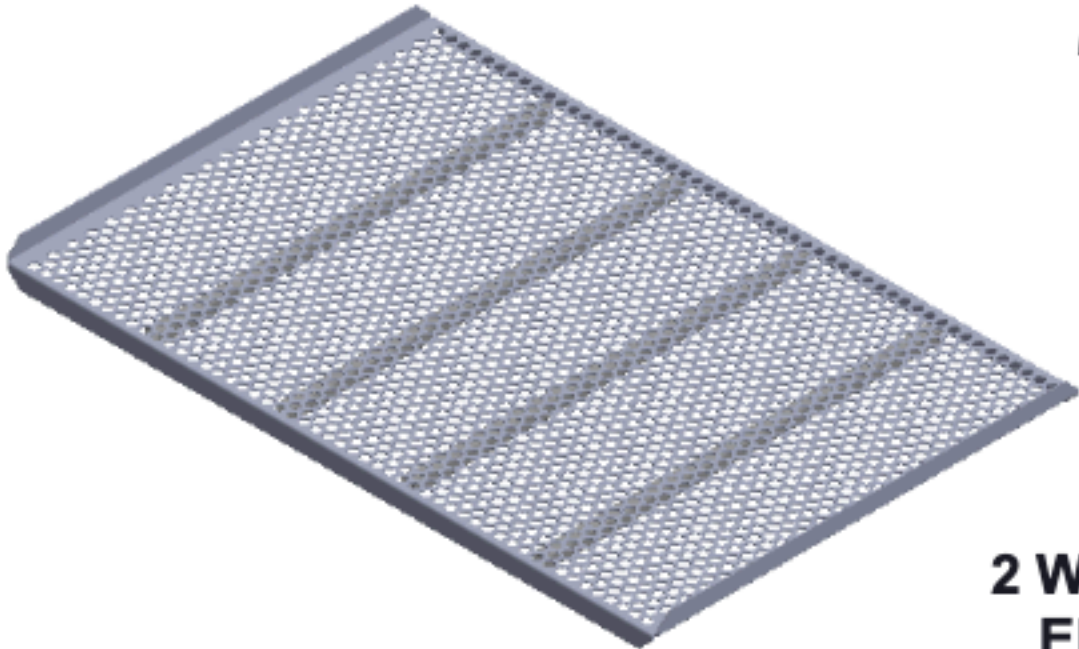


**NOTE: Weld the support tubes to the rack from the top.
Weld the rack's hole edges to the channel flanges.
Grind flush if needed.**

REVISIONS				
REV	CD	DESCRIPTION	ENGINEERING	DATE
-	CO-XXXXXXX	INITIAL RELEASE	DRAWN: F. HUME	02/10/26
			APPROVED:	XX/XX/XX

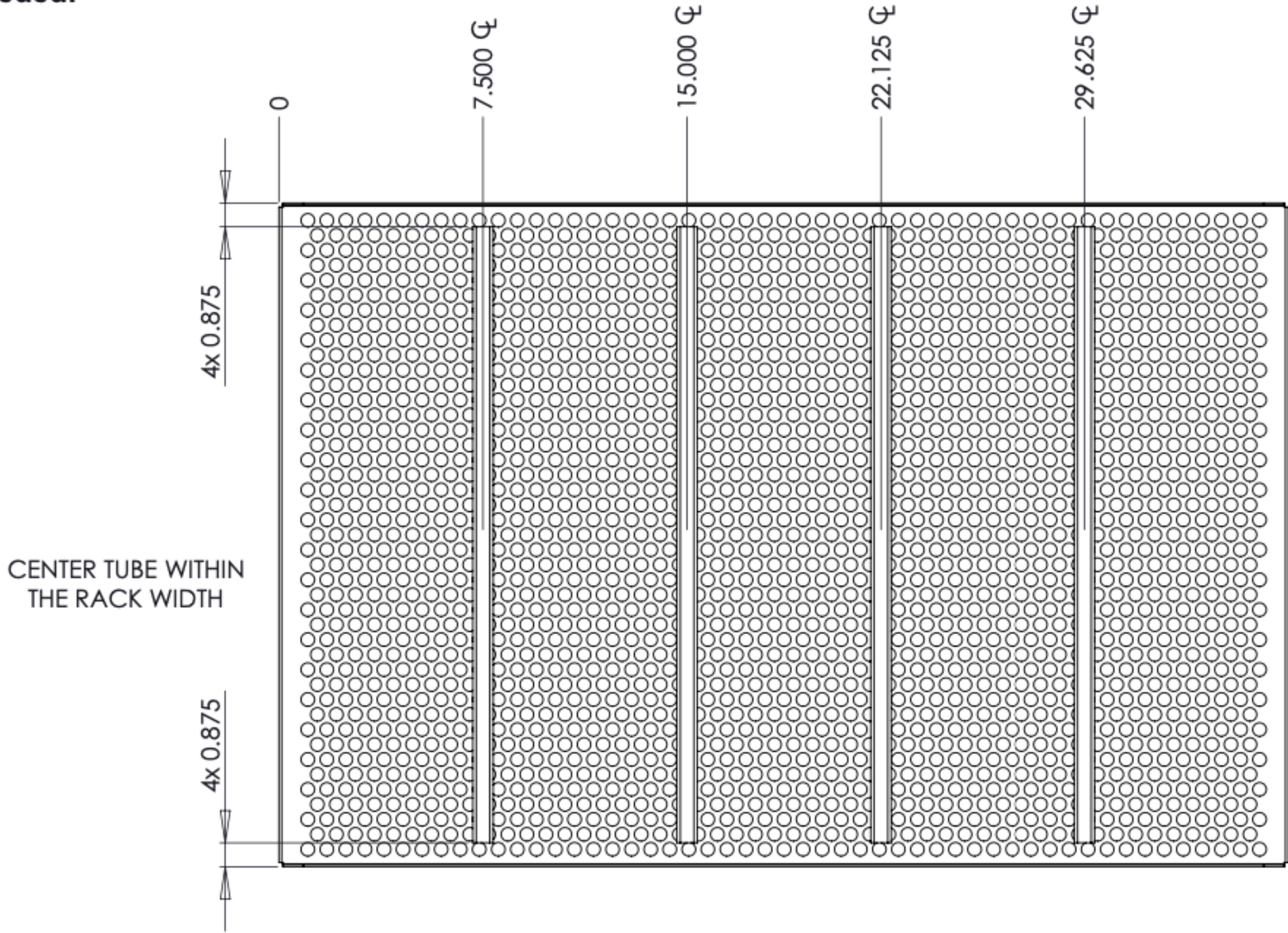

**2 WELDED ASSEMBLIES
ELECTRO POLISHED**

PROPRIETARY INFORMATION THIS DOCUMENT AND THE CONTENTS THEREOF ARE CONSIDERED PROPRIETARY INFORMATION OF ADVANCED STERILIZATION PRODUCTS AND DISCLOSURE TO INDIVIDUALS OR DISSEMINATION, PUBLICATION, OR COPYING IS PROHIBITED WITHOUT PRIOR WRITTEN CONSENT BY ADVANCED STERILIZATION PRODUCTS, IRVINE, CALIFORNIA

**MATERIAL
316 STAINLESS
STEEL**

**FINISH
ELECTRO POLISHED**

INCH

UNLESS OTHERWISE SPECIFIED
DIMENSIONS ARE IN INCHES

TOLERANCES

FRANCTION ± 1/64
ANGULAR ± 30°
O. XX ± 0. 03
O. XXX ± 0. 010

CHAMFER ± 5°
SURFACE FINISH 63 RMS
DEBURR AND BREAK ALL SHARP
EDGES

ASPTM

Bewbodfe Tufsjmj(bujpo Qspevdut

Jswjof- Dbmjgpsojb :3729

**TITLE
Autoclave Rack Assembly, Top & Middle**



B

SCALE: 1:5

SHEET 1 OF 1

REV.
A